

MIRI® GA

Gas and Temperature Validation Unit

MIRI® GA is a tabletop device that is designed to perform easy and continuous external gas and temperature validation of IVF incubators. This device can be used to validate different varieties of Esco Medical IVF incubators, including MIRI® TL, MIRI®, MIRI® II-12 and Mini MIRI® Dry multiroom IVF incubators.



Key Features



Constantly validate up to 6 (six) chambers at the same time.




Controllable flow rate.



Monitor up to 6 x PT1000 sensors.



6 ports for sequential gas samples.

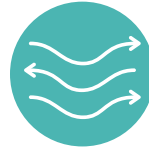

Gas feedback returns sampled gas to incubator or exhaust.



Audible and visual alarms.



As validation unit, MIRI® GA can provide:

- Parameter data logging.
- Adjustable alarm thresholds.

Connected to a PC, the logging software will document parameters and optionally send e-mail alarms if set up.



Want to know more about this product? Scan the QR Code

Esco Micro Pte. Ltd.

19 Changi South Street 1 • Singapore 486 779 • Tel +65 6542 0833
medical@escolifesciences.com • esco-medical.com

Esco Global Offices: Bangladesh | China | Denmark | Germany | Hong Kong | Indonesia | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | UK | USA | Vietnam